

IC30 ACTIVEX™ MATTRESS

The IC30 is our all in one mattress. Incorporating pressure care, bed mobility and most importantly a very comfortable sleeping surface. Highly effective in the treatment and prevention of pressure sores.

RISK CATEGORY:

HIGH RISK

COMFORT CATEGORY:

MEDIUM FEEL



Aged Care




Every Day Care



Disability




Facilities

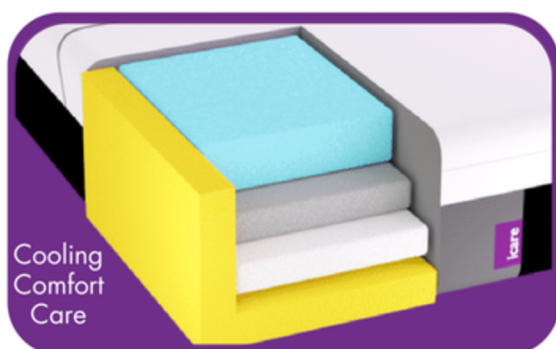

TRANSFERS

Transferring to and from the IC30 mattress is assisted with the XP Support mattress edge and anti-hear BlueTec mattress surface cover.

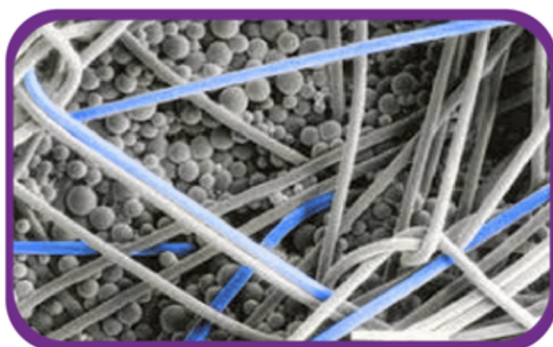

SLEEP QUALITY

Quality sleep is an important factor to wellbeing and healing. The IC30 mattress creates an ideal balance for medical and pressure care needs, sleep quality and longevity.



PRESSURE CARE

The IC30 Mattress is designed to prevent and reduce the risk of pressure injuries. The ActiveX™ sleep surface provides partial immersion with envelopment, creating an excellent average peak pressure, while maintaining bed mobility and lateral movement. The mattress microclimate is regulated by using BlueTec, gel infusion and charcoal.

FEATURES



Pressure Relieving



ZeroTec covers
available



ActiveX™ Gel
temperature regulating



Anti-shear materials



Moisture wicking &
breathable fabric.



ActiveX™ Charcoal
reduces odors and
is microbial.



Multi layer
immersion



Reduces Rollout
Risk




Contour-Latch™
attaches to most
profiling beds.




20+ waterlow
rating

BED CODES & SIZES

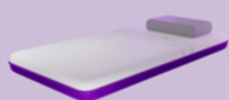


SINGLE
Code: IC30S
Mattress Size: 1900 x 900

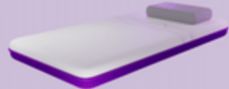


LONG SINGLE
Code: IC30LS
Mattress Size: 2030 x 900




KING SINGLE
Code: IC30KS
Mattress Size: 2030 x 1070

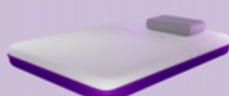


DOUBLE
Code: IC30D
Mattress Size: 1900 x 1350

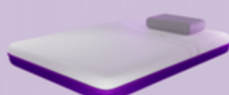

LONG DOUBLE
Code: IC30LD
Mattress Size: 2030 x 1350

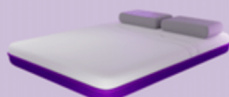

QUEEN
Code: IC30Q
Mattress Size: 2030 x 1520

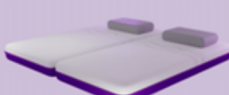

SPLIT QUEEN
Code: IC30SQ
Mattress Size: 2030 x 750 (per side)

SPECIFICATIONS

ActiveX™ pressure relieving material
Vant ventilated sides and base fabric
Bluetec temperature regulating cover
Option Zerotec waterproof cover (see pg. 22)
Flip free design
Zip off machine washable covers
30kg - 200kg user weight
Standard cover MVTR: NA.
Medical Cover MVTR (g/24h-m2): 1635.0 Average
Medical Cover BS7175: Crib 5
Waterlow ≤20

Feel Rating	Medium
Mattress Height	200mm
Risk Category	High Risk
User Weight	20kg - 228kg
Warranty	10 Year Warranty

PRESSURE CARE

Risk Category High

This rating is based on our average peak pressure of <19mmHg. For some high risk clients Icare Zerotec covers are recommended.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For some high risk clients Icare Zerotec covers are recommended. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.

CARE & CLEANING

To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

NEW


IC30 ACTIVEX™

icare
Medical Group

The IC30 Gel-Infused ActiveX™ Mattress with Charcoal and Firm Edge Support is designed to provide comfort and support for a good night's sleep. The gel infusion in the mattress helps to regulate temperature, ensuring a cool and comfortable sleep surface. The addition of charcoal helps to absorb odors and moisture, promoting a clean and hygienic sleeping environment.

ActiveX™ is temperature and pressure sensitive which provides an optimum balance of immersion and envelopment creating a sleep surface with excellent support, comfort and very high level of pressure care.

The firm edge support helps to maintain the shape and structure of the mattress, reducing the risk of roll-out and assisting with transfer access. This also makes it ideal for use on profiling beds, which are beds that can be adjusted to various positions to provide comfort and support for those with limited mobility.

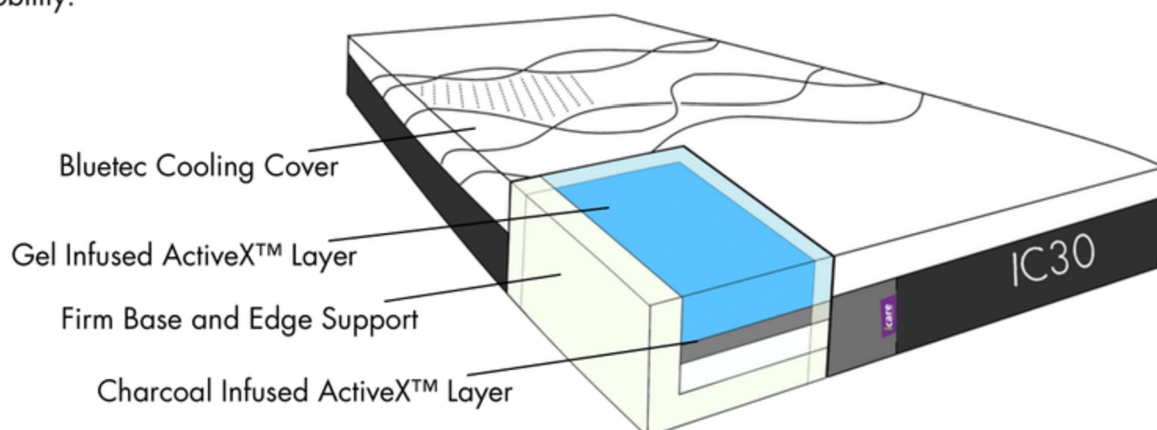

SPECIFICATIONS

Feel Rating	Medium
Mattress Height	200mm
Risk Category	High Risk
User Weight	20kg - 228kg
Warranty	10 Year Warranty

BED CODES & SIZES

LONG SINGLE	IC30LS	2030 X 900mm
KING SINGLE	IC30KS	2030 x 1070mm
LONG DOUBLE	IC30LD	2030 x 1350mm
QUEEN	IC30Q	2030 x 1520mm

PRESSURE CARE

Risk Category High

This rating is based on our average peak pressure of <16.5mmHg.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For more information refer to the Icare pressure care manual.



movemobility.com.au

IC30 ACTIVEX™ MATTRESS



Pressure Redistribution Information

IC30 MATTRESS


Risk Level: High
SWL Minimum 20kg
SWL Maximum 228kg

Waterlow ≤ 20
Immersion Rating: Medium
Envelopment Score: 4/5


Subject 1: Female. 58kg
Average Peak Pressure: 14.3mmHg
Support area: 2379cm²



Subject 2: Male. 84kg
Average Peak Pressure: 18.7mmHg
Support area: 2548cm²

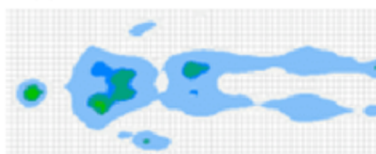

Scale mmHg

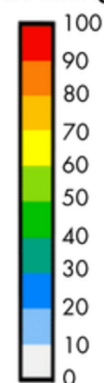

Risk category ratings, immersion ratings and envelopment scores are based on average peak pressure, therapist recommendations, case studies and Waterlow scale. For some high risk clients Icare Zerotec covers are recommended. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of five days.