



novis  **ORNAMIN**

Ornamin

Smart Design brings everyone to the Table



MADE IN GERMANY

A RANGE OF ERGONOMICALLY DESIGNED EATING AND DRINKING AIDS FROM ORNAMIN

Developed in consultation with clinical experts and customer needs, this comprehensive range of eating and drinking aids is both functional and elegant.

Specifically designed for people suffering from Alzheimer's (a form of dementia), this colourful range of tableware is ideal for anyone with cognitive impairment or decreased visual capacity. Many items are also useful for people with impaired mobility or reduced motor skills such as Parkinson's Disease.

In addition to symptoms such as forgetfulness, speech impediments and orientation difficulties, people with Alzheimer's struggle to distinguish items through colour contrast. Along with other issues such as depression or restlessness, this can result in a decreased appetite and a reduced consumption of both food and drink, which can lead to significant weight loss and dehydration.

Studies in clinical nutrition have shown a significant increase in both food and liquid intake when using tableware that offers a high colour contrast (when compared to the table or the food) – with red offering the highest contrast and therefore the most successful result.

Ornamin have applied these results to develop a comprehensive and functional range of contrasting eating and drinking aids, designed to help stimulate the appetite and orientate anyone impacted by contrast sensitivities. Made from lightweight yet sturdy melamine, the colour contrasting designs are easily recognisable and perfect for independent dining.



ORNAMIN
family tableware




VITAL SERIES

A range of red and white contrasting tableware designed to stimulate the appetite and help people more easily see their food or drink due to the colour contrasting rims.

Each piece is designed for maximum functionality and ease of use, making it easy for anyone with cognitive impairment or other mobility limitations to dine with dignity and independence.

Vital Plate and Bowl

- ❑ sloping inner edge helps food be easily scooped
- ❑ inconspicuous overhang for added grip and stability
- ❑ anti-slip ring on base to stop unwanted movement
- ❑ ideal for single handed use
- ❑ plate edge higher than a standard plate
- ❑ plates suited to main meals and soups
- ❑ bowl suited to entrees and desserts, or for smaller portions

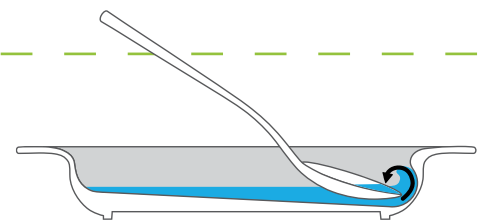


CROSS SECTION

Vital Ergo Mug

- ❑ conical shaped inner cup
- ❑ two broad handles for easy grip and added stability
- ❑ shape allows contents to be emptied without tilting head
- ❑ ideal for drinking in bed or laying down
- ❑ thermodynamic function helps contents stay warm or cool

IMPEDIMENTS WITH CONVENTIONAL MUGS

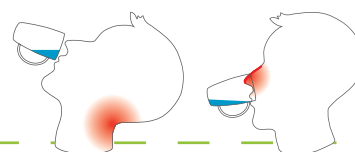

VITAL ERGO MUG

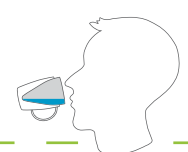


PLATE 26 CM



PLATE 20 CM



BOWL 15 CM



ERGO MUG

	PLATE 26 CM	PLATE 20 CM	BOWL 15 CM	ERGO MUG
CODE	10521	10144	9657	12489
CAPACITY	760 ml	360 ml	185 ml	200 ml
WEIGHT	495 g	300 g	155 g	147 g
DIMENSIONS	260 x 40 mm	200 x 15 mm	155 x 35 mm	80.5 x 108.5 mm



THERMO SERIES

An innovative range of tableware as simple as it is ingenious. Designed to keep food warmer for longer, these products are ideal for slow eaters, ensuring that the meal or drink stays warm and flavoursome right to the last mouthful.

Through a special opening in the outer wall, the plate or bowl can be filled with hot water to help keep the food warm. At the same time, it helps reduce the risk of burns while keeping the heat on the inside where it is needed. Alternatively, the tableware can also be filled with chilled water to help keep food items cold as needed.

Made from a high-quality melamine inner with an external thermoplastic material outer and two grip inserts made of a soft, non-slip silicone plastic. The handles provide a firm hold while also covering the fill holes to prevent the hot or cold water from spilling.

Thermo Plate and Bowl

- ❑ thermodynamic function helps contents stay warm or cool
- ❑ red and white colour for maximum contrast
- ❑ extra deep bowl and broad plate
- ❑ non-slip handles for added grip and stability
- ❑ anti-slip ring on base to stop unwanted movement
- ❑ replacement handles available (sold as pair)



CROSS SECTION: FILL THE BASE WITH HOT OR COLD



GLOW IN THE DARK



CROSS SECTION MUG



CROSS SECTION MUG VITAL

Thermo Mugs

- ❑ thermodynamic function helps contents stay warm or cool
- ❑ large ergonomic handle for easy grip
- ❑ white interior ensures colour of drink is maintained
- ❑ Vital Thermo Mug offers conical shaped inner cup
- ❑ Glow-in-the-dark Mug glows after exposure to light



PLATE



BOWL



MUG



MUG VITAL



GLOW



HANDLES

	PLATE	BOWL	MUG	MUG VITAL	GLOW	HANDLES
CODE	10863	10858	11797	11792	5124	9726
CAPACITY	425 ml	300 ml	320 ml	200 ml	280 ml	-
WEIGHT	450 g	224 g	155 g	145 g	115 g	15 g
DIMENSIONS	255 x 42 mm	190 x 63 mm	80 x 108 mm	80 x 108 mm	71 x 115 mm	16 x 85 mm

ERGO PLUS CUP



Ideal for people with tremors or limited hand dexterity, this sturdy cup has two wide handles for safe handling and stability. Handles can be grasped by the hand or slid over the hand for added security when hands are placed around the cup.

Available in a range of colours, the Ergo Plus includes a measurement scale on the inside of the cup to help with easy volume control.

Use with the Non-Drip Spouted Lid for added safety against spillage.

■ Weight 90 g ■ Capacity 250 ml ■ Dimensions 71 x 100 mm



COLOUR	NATURAL	YELLOW	ORANGE	RED	BLUE
CODE	3187	3190	3191	3194	3195

NON-SLIP CUP



More than just a colourful design, these innovative Non-Slip Cups are made from a special non-slip material to ensure a firm grip even when hands are damp or trembling. Available in a range of colour combinations, the contrasting designs help to orientate the hand to the middle of the cup.

The Non-Slip Cup comes in two bright designs – a Flower or Measuring Icon design. The vibrant colours of the Flower design make it especially attractive to people suffering from visual or cognitive impairment, while the Measuring Icons make it especially easy to pour and record volume measurements. The measurement amount (in millilitres) is highlighted in the centre of each circle.

Use with the Non-Drip Spouted Lid for added safety against spillage.

■ Weight 42 g ■ Capacity 250 ml ■ Dimensions 71 x 100 mm



	TRANSLUCENT FLOWER				OPAQUE FLOWER				MEASUREMENT			
COLOUR	WHITE	ORANGE	GREEN	BLUE	BLUE	YELLOW	GREEN	RED	BLUE	GREEN	ORANGE	WHITE
CODE	8713	8707	8709	8711	9484	9486	9487	9488	5941	5943	5945	5946

NON-DRIP SPOUDED LID



An optional accessory to suit most mugs or cups in the Ornamin range, this spouted lid can be sealed over the top of the cup to prevent leaks or spills.

The spout also makes it easier to drink without the need to tilt the head, and helps to control the flow of liquid. The user is required to gently suck to release the fluid, so controlling the rate of fluid release and minimising the risk of choking.

Promotes independence for anyone with limited mobility or dexterity, or anyone who struggles with ordinary drinking function following a stroke or other debilitating illness.

CODE 6972 ■ Weight 20 g

NON-SLIP PLACEMATS

Made of durable polypropylene, these non-slip placemats have a special coating on both sides of the mat to help when eating and ensure against unwanted slippage.

The anti-slip property of the top surface is guaranteed to hold dishes firmly and securely, ensuring the plate does not need to be held while eating. The special coating on the bottom surface prevents the mat from slipping on the table or eating surface. Easy to wipe clean.

□ Dimensions 395 x 275 mm



COLOUR	BLUE	RED	YELLOW	GREEN
CODE	6206	6208	6209	6207

NON-SLIP BOARD

The simple yet innovative design of this melamine non-slip board allows people with limited mobility or hand dexterity to prepare food or butter bread without the need for help.


The clever design suits both right and left handed users.

Three sides of the board have a raised edge to hold bread or other food in place and prevent it from moving. The base on the board has an anti-slip ring to ensure the board remains secure on the preparation surface. These features make it easy to butter bread with just one hand.

Available in a range of bright colours, these boards can be stacked for optimum space saving in the kitchen. Dishwasher safe.

Use with the optional Food Preparation Aid to help hold fruit, vegetable or cheese.

□ Weight 360 g □ Dimensions 220 x 170 mm

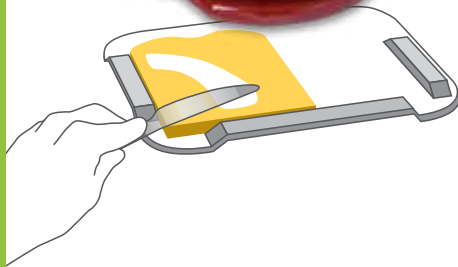


COLOUR	BLUE	RED	YELLOW	GREEN
CODE	6599	6597	6598	9109

FOOD PREPARATION AID

An optional accessory for the Non-Slip Board, this useful tool can be used to help when cutting food such as fruit, vegetables or cheese. Ideal for people with limited mobility or dexterity, or only the use of one hand.

The small plate with four spikes attaches easily to the edge of the Non-Slip board, and helps to hold whatever small item is being cut, peeled or prepared.

CODE 11449 □ Weight 35 g

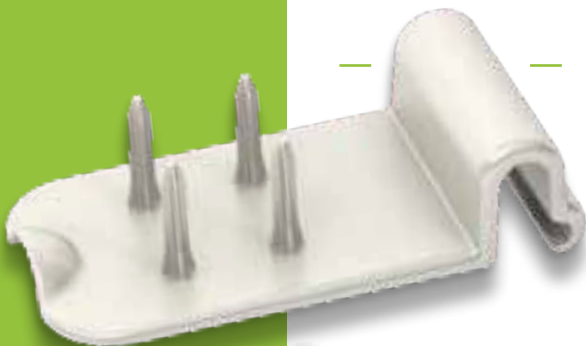

PRODUCT FINDER

FINDING THE BEST SUPPORT – THE APPROPRIATE HELP FOR EVERY NEED

The various eating and drinking aids by ORNAMIN provide targeted aid for specific needs. This overview shows you which products offer support for which symptom.

SYMPTOM	ITEM	
	Non-Slip Board	
Limited hand movement or shaky hand movements	Plate Vital Ø 26 cm	
Impaired motor skills	Plate Vital Ø 20 cm	
Poor or limited vision	Bowl Vital Ø 15 cm	
	Non-Slip Placemat	
Slow eaters	Thermoplate Thermobowl	
Restricted of movement or paralysis of the upper extremities and the neck area	Vital Ergo Mug Thermo Mug Vital	
Slow eaters and drinkers	Thermo Mug Thermo Mug Vital	
Numbness in the fingers, moist or shaky hands	Non-Slip Cup	
Shaky hands such as parkinson's disease and dementia	Ergo Plus Cup	
Difficulty swallowing, shaky hands, poor or limited vision	Non-Drip Spouted Lid	
Dangers of choking, slight limitations in the neck		



Eating together is an important social event and one of life's simple daily pleasures – provided everyone can participate and enjoy their food. The award-winning family tableware by ORNAMIN has been designed with style as well as function, removing any stigma or embarrassment about using adaptive crockery because its supportive features are hidden in its design.

novis


ORNAMIN

1300 738 885 / novis@novis.com.au
12/12 Mars Rd, Lane Cove West, NSW 2066

Scan QR code for more information
novis.com.au/ornamin



08042021F